FX4000

MANUALE D'ISTRUZIONE

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FX4000

INSTRUCTION MANUAL

Congratulations! You have chosen a quality thermostimulator designed with utmost care, among the most efficient available on today's market. The careful design, its full compliance with current safety standards and its ease of use place it at the forefront of the equipment currently manufactured in this sector.

GENERAL WARNINGS

- To identify the instructions for the model of thermostimulator in your possession, read the data on the label affixed to the rear.
- Before connecting up the apparatus to the main supply, check that the data indicated on the label (voltage and absorbed power) correspond to the characteristics of the electric power system.
- Before mounting the apparatus on the stand or the wall-mounted arm, be sure it is disconnected from the electric power supply
- Make sure not to rest towels, cloths, etc. on top of the thermostimulator, but leave the apparatus free so that air can circulate through the slits of the grill provided on the top; without proper air intake, the apparatus will be unable to work correctly.
- The Company declines all responsibility for any injury to persons or damage to things caused by installation or assembly that not been carried out in compliance with the General Warnings, or using supports and/or accessories other than those manufactured by us, or for any injury to person or damage to things caused by installation or assembly that have been carried out by non-specialized or non-authorized staff.
- In the case of the wall-mounted arm, the Company moreover declines all responsibility for injury or damage resulting from the unsuitability of any support made of masonry or other material.
- The apparatus can be cleaned using a cloth moistened with water or ordinary commercially available detergents, provided they do not contain alcohol or acetone. However before cleaning, make sure that the appliance is disconnected from the electric power supply and that it has cooled down properly.
- For replacement of internal parts, as well as in the event of damage to the power supply cable. In all cases make sure to call in specialized and authorized staff.
- If the power supply cable is damaged, it must be replaced from a cable or a together special available near the constructor or its service technical attendance.

ATTENTION

THE HEATING ELEMENTS (LAMPDS) DO NOT HAVE THEMSELVES TO BE ILLUMINATED DURING THE OPERATION

TREATMENT TABLE

The times, temperatures and ventilation recommended are purely indicative. They may vary according to the type of hair to be treated, the product used and the temperature set. Give the wide range of products and all the different kinds of hair the hair stylist must be able to assess correctly the time required and the temperature that is most suitable for each individual treatment. Consequently, the manufacturer declines all responsibility for any damage that may derive from improper use of the equipment.

ATTENTION - Never use this appliance near water container in bath tubs, sinks or other containers. - This hair dryer must be fixed far from people using the bath tub or the shower.



- Protect the eyes from cancellations IR, don't fix the lamps directly if these are ignited. In case of extended exposure of the eyes to the direct light of the lamps, consult a doctor.

IMPORTANT - water poses a serious hazard even when the devices is switched off; do not leave device connected to power supply after use in the bathroom. - For additional electrical protection the installation of a differential power unit (max. 30mA) is recommended. - For replacement of internal parts or in case of damage to the power cord, contact the manufacturer or service centre, or similarly authorized personnel. To the aim to assure an additional protection it is advised to install, in the electrical circuit, who feeds the premises from bath, a differential

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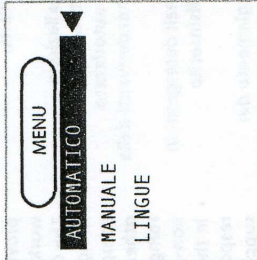
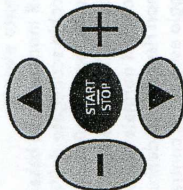
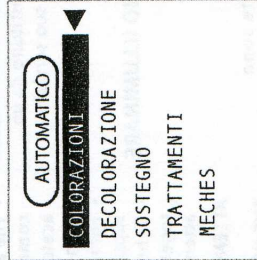
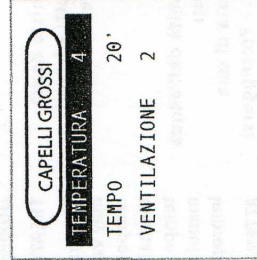
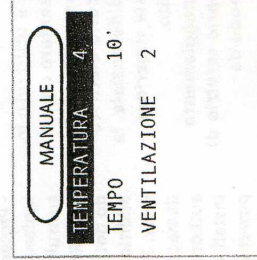
The device is switched on by pressing the button indicated with **O / I** (A). At the screen containing the name and corporate logo and name of the product.


Choice of language. Will I choose the desired language pressing the selection keys. Push button more (C) to **MENU**.

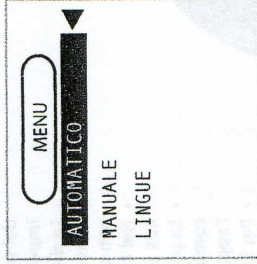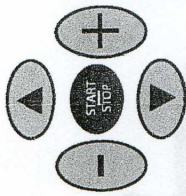
MENU screen. Will I choose between three options going down or up with the keys selection (B - D). You can go to **AUTOMATICO** menu pressing more key (C).

Going over the other two options and pressing button more (C) you can access to relative screens.

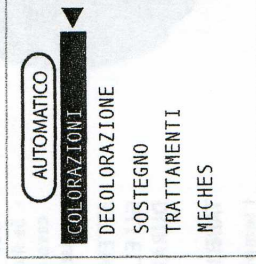
AUTOMATICO screen. Will I choose between five options going down or going up with the keys selection (B - D). Access to the screen or other desired pressing the button more (C).

Will I return to the screen **MENU** going to **AUTOMATICO** and pressing less key (F). Going on the other two options and pressing key more you will access to their screens.

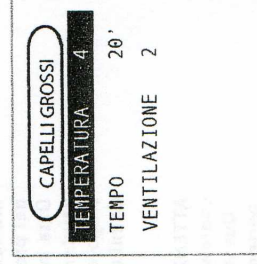Screen **PROGRAM SELECTED**. Going on **TEMPERATURE**, **TIME** OR **VENTILATION** with the directional buttons (D-B) you can change the numerical values using the more or less keys (F-C), or can leave the set values. P.S. editing a numerical value the current screen pass to **MANUAL**. Pressing **START / STOP (E)** the program will start. Going on the **name of program** with directional arrows and pressign the less key (F), will return to the previous screen.

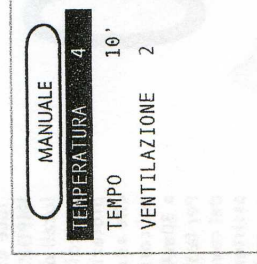Screen **MANUAL**. On **TEMPERATURE**, **TIME** or **VENTILATION**, you can change the values using the more and less keys (C-F). The **TEMPERATURE** will have values ranging from 1 to 6. The **TIME** will have values ranging from 1 to 60 minutes. The **VENTILATION** have values ranging from 1 to 3. The program will be lstarted pressing **START / STOP** key(E). Going on **MANUAL** with the directional buttons (D-B) and pressing the key less (F) you'll


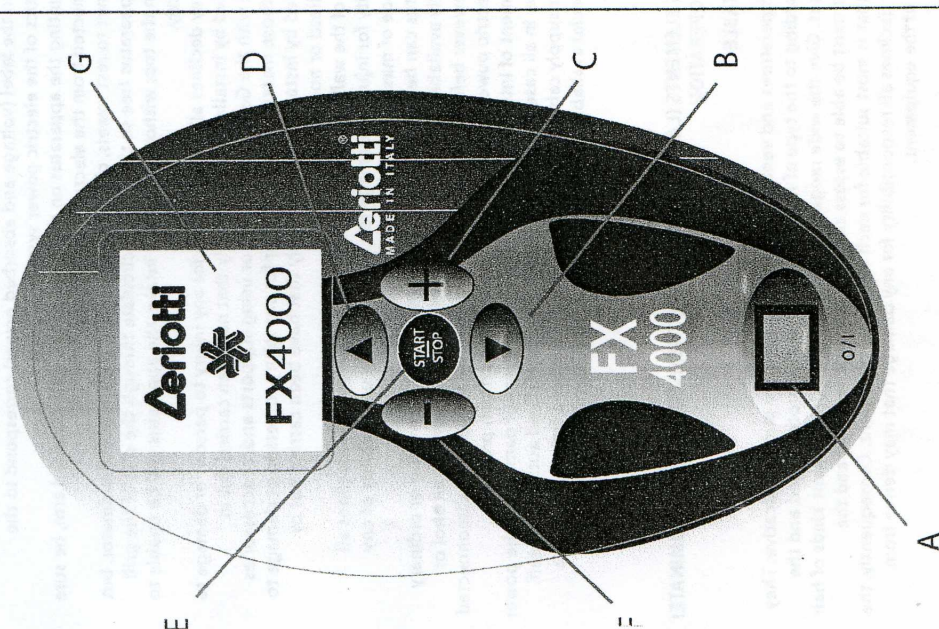
Fig. 1

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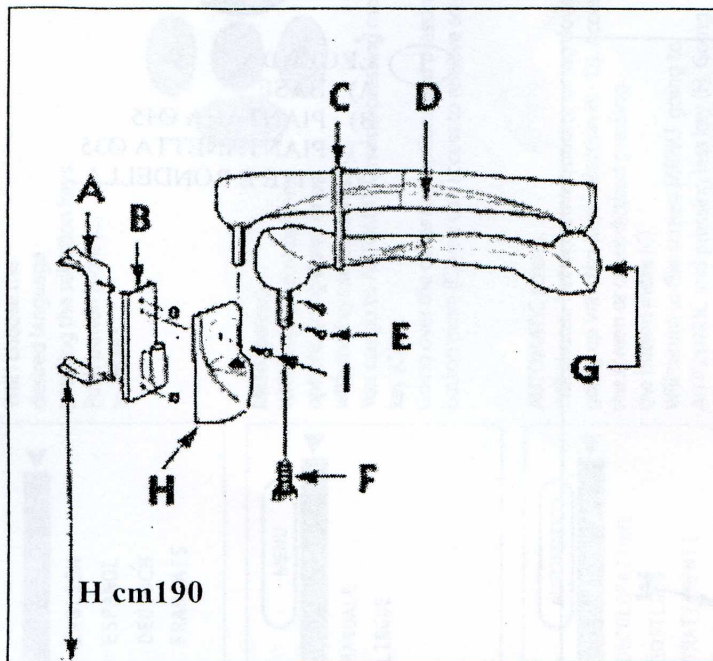
GB WALLARM – INSTRUCTION FOR ASSEMBLY

1. Set the bracket (A) in the wall; recommended height 190 cm. Max.
2. Fasten the plate (B) to the bracket (A) using the two nuts.
3. Mount the pact cover (H) on the plate (B) and fasten with the screw (I).
4. **IMPORTANT:** Without removing the safety clamp (C), witch withholds the articulated part of the arm inside which the return spring is bossed, insert the arm (D) on the plate (B)
5. After loosening the two screws (E) at the end f the arm (D), insert the pin (F) of the fork and tighten the screws right home.
6. Only now's it possible to remove the safety clamp, holding the apparatus with one hand to prevent any Kick-back in the working position.
7. Connect the plug of the apparatus to the socket provided on the arm and connect up to the electric power supply. Should use cause an alteration in the balancing mechanism, correct it as follows: if lifting the apparatus the counterweight tends to drop, turn clockwise the screw (G) on the elbow of the arm clockwise, using a n.8 Allen wrench; if the hood tends to do the opposite, turn it counter clockwise.

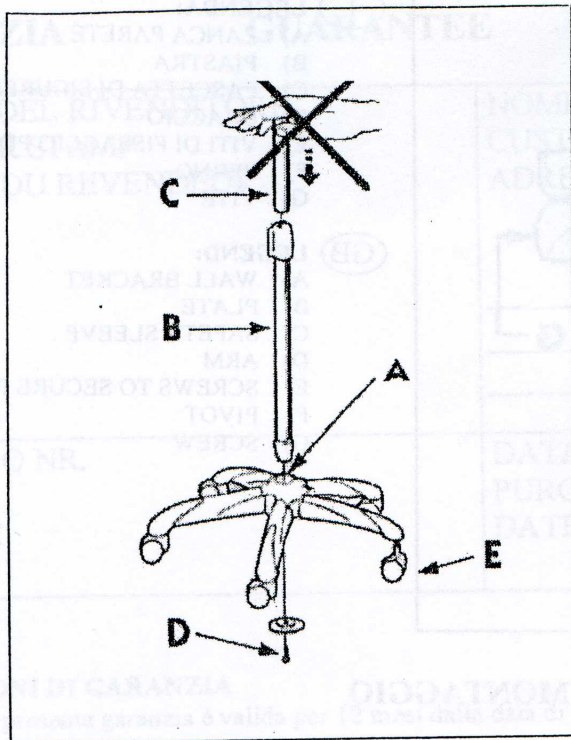- I LEGENDA:**
- A) BASE
 - B) PIANTANA Ø45
 - C) PIANTANETTA Ø35
 - D) VITE E RONDELLA

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
GB FLOOR BASE - ASSEMBLY INSTRUCTIONS

1. Introduce the supporting pole (B) into its housing at the top of the base (A), then secure it by tightening the screw with washer (D) at the bottom.
2. Loosen the plastic sleeve of the pedestal (B), turning it counter clockwise. Insert the small pedestal (C), resting it in the socket. Make absolutely sure not to compress the spring inside the pedestal (B): if this is pressed too hard, it could push out the small pedestal, with possible serious consequences. Fasten the small pedestal (C) on, turning the sleeve clockwise.
3. Insert the appliance onto the stud (C) and regulate its height, turning the sleeve as described in point 2.

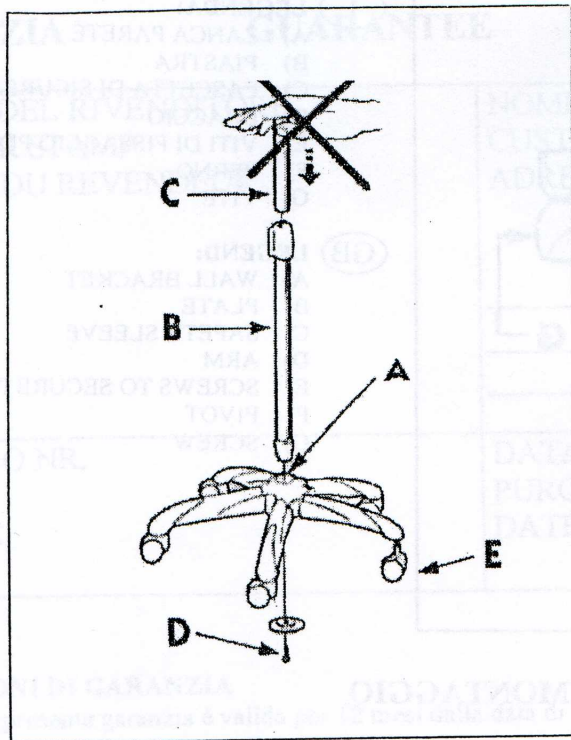- I LEGENDA:**
- A) BASE
 - B) PIANTANA Ø45
 - C) PIANTANETTA Ø35
 - D) VITE E RONDELLA

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WARRANTY CONDITIONS

- 1) This warranty is valid for 12 months from the date of purchase, reported on this certificate and on the fiscal receipt.
- 2) Within such period the Company undertakes to repair or substitute, without charge, any part showing manufacturing or material flaws
- 3) This warranty does not cover the following: **A-** Damage caused by incorrect use or maintenance, or by normal wear. **B-** Malfunctions caused by unauthorized personnel tampering. **C-** Damage caused by mishandling, accident or transport.
- 4) In the period of validity of this warranty, all covered repair and substitutions shall be performed exclusively by the manufacturer. All shipping expenses shall be charged to the buyer.
- 5) Whenever you request the application of this warranty, enclose this certificate, duly filled, and a copy of the fiscal receipt, with the appliance.
- 6) Conserve the packing case for some days after purchase

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